

Pura



Multifunctional chair

Variable solutions for One Day



Pura Day Surgery

Fully equipped chair for day surgery, recovery or diagnostic procedures.

Care



Pura Treatment

Comfortable for dialysis, oncology treatment or blood donation.

360° care at Pura chair

TRANSPORT IS EASY

Transporting the patient is safe and is not physically strenuous.



PLEASANT STAY

The comfort of the patient is ensured during long-term procedures by ergonomical frame and soft padding.



FOCUSED ON SAFETY

Unique siderails design ensures patients safety while Soft Drop mechanism enables for smooth & easy lowering of siderails and easy access to the patient.



RECOVER AND GO

The low, safe position is comfortable when staying in or when entering or exiting the chair.

Focused on safety

The new siderail concept of the Pura Surgery chair increases the safety level of professional outpatient care and provides patients with increased protection and comfort. The applied LINET technologies support ergonomics and greatly facilitate the work of nurses.

Soft Drop technology

The siderails are equipped with Soft Drop technology for optimising and softening the movement. The smooth and quiet lowering of the siderails is comfortable, not only for the staff, but it also makes the patient's stay in the chair more pleasant.




Trendelenburg position

The anti-shock position, or the so-called Trendelenburg position, is used for short-term collapse and is one of the pre-set functions on the Pura chair.



Mechanical CPR

A quick mechanism allows for the simple adjustment of optimal CPR conditions.

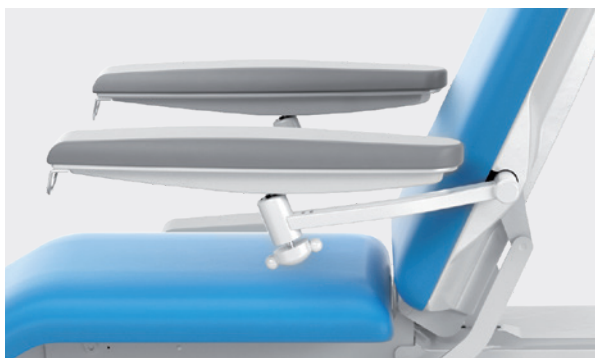

Ergonomical armrest

The Pura Treatment chair is equipped with comfortable armrests for a long-term stay during dialysis or oncology infusion therapy.

Variability in positioning

The quick and easy electric control adjusts the height of the chair for ergonomic work with the patient.

The variable range of positions of the back and leg section allows the Pura chair to be used in a wide variety of ways, from therapeutic procedures in the horizontal position to individualised adjustment and easy adjustments during long stays in the chair for treatment procedures.



Fowler's position

The Pura chair allows you to use the so-called Fowler's position to optimally place the patient in a semi-recumbent position to relieve pain and support breathing.



Portable handset

The user-friendly and ergonomic handset includes a pre-set TR, flat and welcome position. Also a safety lock of the functions, GO button and other intuitive features for the easy and quick adjustment of the chair.

Transport is easy

The driving and manoeuvrability features allow for fast and safe patient transport with minimal physical exertion.

Easy and safe transport

The state-of-the-art structure of the undercarriage fitted with wheel with a 150 cm diameter castor and a central brake is the best option for transporting patients in a chair. Along with its well-designed details, such as the transport handle or anti-crash bumpers, it creates excellent conditions for driving in outpatient facilities and hospitals.



Transport handle

The foot section can be equipped with a robust handle that allows ergonomic handling of the chair, especially during transport.



5th castor

A 5th castor on the undercarriage is an optional feature. This ingenious mechanism allows you to better manoeuvre the chair in confined spaces, corridors, or in lifts.

Safe and comfortable care

To maximally facilitate the work of nurses and staff, the Pura chair introduces an innovative infection control concept and a sophisticated set of practical details.



Infection control

The clean column construction line together with the compact surfaces form the basis for easy hygiene.




Oxygen bottle holder

It is designed for holding the oxygen bottle in a horizontal position.



Plastic box

Space for personal things.

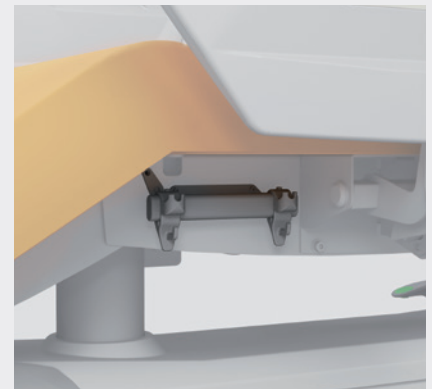

Urine bag holder

It serves for the lodgment of the urine bag (maximally 2 l).

The design of the chair has been primarily adapted to make sure that the patient's stay is pleasant. Soft padding and the material used for the upholstery are the basic prerequisite.



Entering and exiting the chair

The low, safe position is comfortable for the patient, both when staying in the chair and when entering or exiting the chair.



Lockable box for personal things

It is designed for locking small patient's personal effects such as keys, telephone, wallet. The box includes keys for the lock.



Paper roll holder

This multifunctional chair is equipped with integrated paper roll holder in head board.



Removable table

Practical table with separate holder increasing patient comfort.

Technical parameters

Pura Treatment

- 1 Armrests
- 2 Footrest
- 3 Infusion holder
- 4 Transport handle



Technical parameters

Length	Total rest area	2,123 mm
Width	Seat section	590 mm
Height (castors 100 mm)	MAX. – rest area to seat section	934 mm
Height (castors 150 mm)	MAX. – rest area to seat section	970 mm
Castors	Castor diameter / brake	100 mm / individually braked
	Castor diameter / brake	150 mm / central brake
Adjustment angles	Back section	70°
	Armrest – vertical	+ 42° / – 8.5°
	Armrest – horizontal	360°
Electrical parameters	Protection against water and dust	IPX4
	Device protection class	I.
	Classification of applied parts	B
	Batteries back up voltage	24 V / 1.2Ah
	Voltage – input	100–240 V
Trendelenburg position (TR)		12°
Safe Working load		205 kg
Maximum patient weight		190 kg

Pura Day Surgery



1 Siderails

3 Infusion holder

2 Castors

4 Transport handle

COLOUR OPTIONS

Standard colours



Blue



Yellow



Orange



Green



Lavender

Non-standard colours



Pura

Patient

- New safe siderails
- Comfortable access and exiting
- Self-sufficiency and comfort



- Higher degree of functionality and better ergonomics
- Above-standard transport possibilities
- Easy cleaning and hygiene



Staff

Management

- A design for higher prestige
- Excellent quality-price ratio



 wissner-
bosserhoff

LINET

Members of LINET Group

LINET spol. s r.o.

Želčovice 5 | 274 01 Slaný | Czech Republic

tel.: +420 312 576 400 | fax: +420 312 522 668 | e-mail: info@linet.com | www.linet.com



Manufacturer: Borcad Medical a. s., 739 45 Fryčovice 673, Czech Republic, Member of LINET Group

LINET | Art. n° 9601-1070-01 | Edition 05/2021 | Different colour shades possible. The right to change technical parameters is reserved.